



CORAIL⁺

BALLOON-TIPPED GUIDE CATHETER 6F AND 8F - BRAIDED

The inner braiding of the CORAIL⁺ catheter gives it a very high resistance to kinking. Available in 6F or 8F, its distal balloon can perform a temporary occlusion while its large inner lumen allows an easy navigation for another catheter.

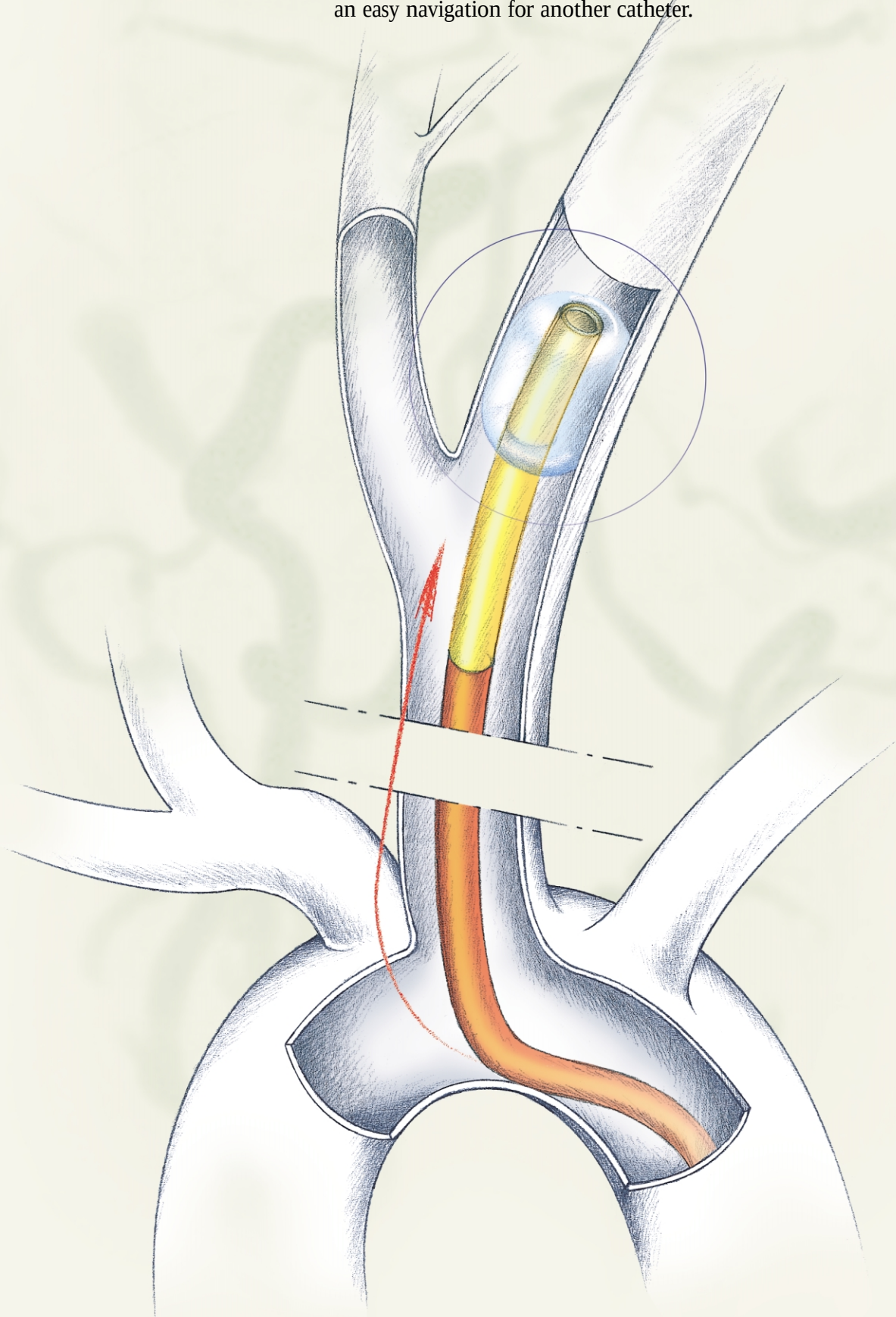


BALT
EXTRUSION

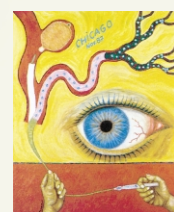

CORAIL+ BRAIDED


CORAIL 6F+

Balloon size Volume (ml)	Ø (mm)
0,1	5
0,2	7
0,3	8
0,5	9
0,7	10
1	11

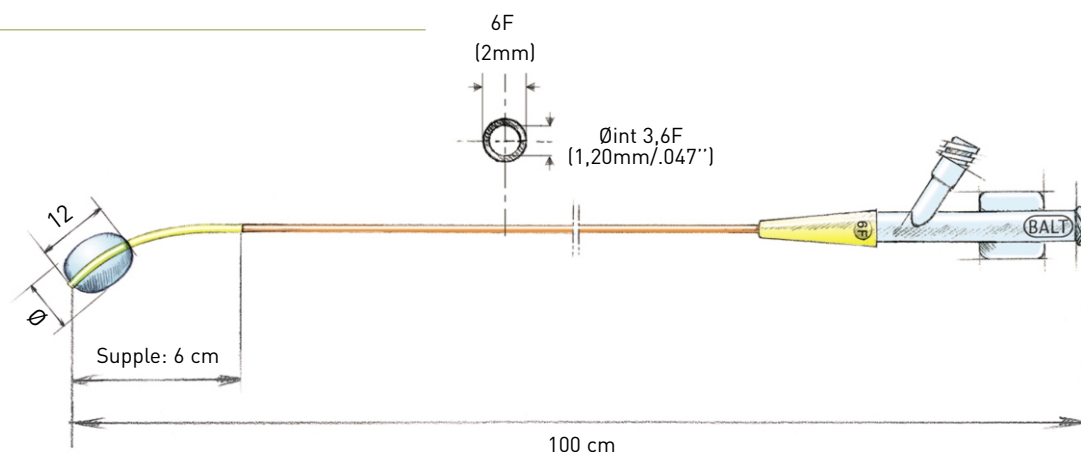

CORAIL 8F+

Balloon size Volume (ml)	Ø (mm)
0,1	6
0,2	7
0,3	8
0,5	9
0,7	10
1	12

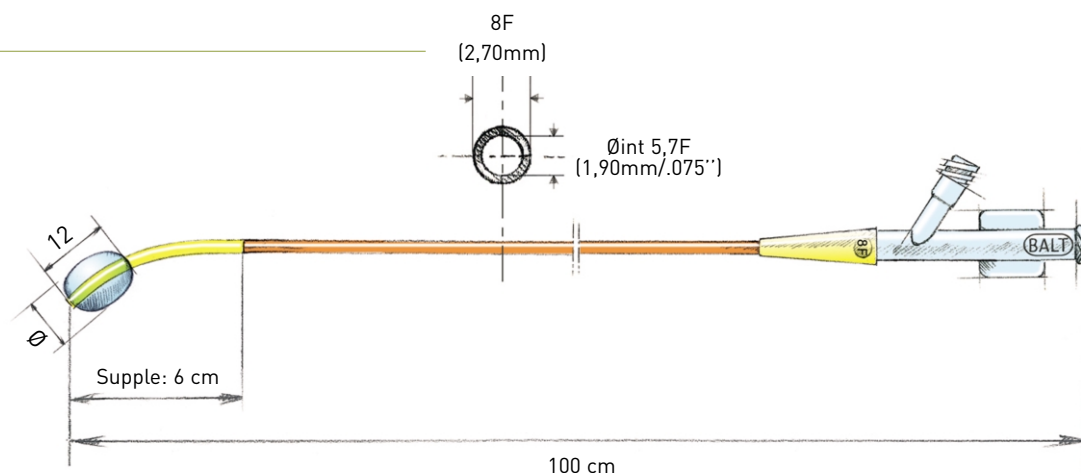

INDICATIONS

- Introduction of detachable balloon (GOLDBAL)
- Clamping test
- Control of the blood flow during embolisation (ethanol...)
- Thrombo-aspiration and/or implantation of stent (Wall stent® type) for acute treatment of carotid dissections

MAIN CHARACTERISTICS:

Reference	O.D.	I.D.	Inflated balloon 1 ml Ø x L	Supple part	Total
CORAIL 6F+	6F (2mm)	3,6F (1,20mm/.047'')	11x12	6 cm	100 cm
CORAIL 8F+	8F (2,8mm)	5,7F (1,90mm/.075'')	12x12	6 cm	100 cm